

Step 1 - Collect Blood Sample




- ✓ Complete Requisition Form
- ✓ Fill tubes to capacity, refer to the tube's reverse side
- ✓ Label vials: Full legal name, DOB, Date & Time of Collection.

Step 2 - Package Correctly




- 1 Place TUBES in SLEEVE > Place in BIOHAZARD BAG
- 2 Place REQUISITION FORM into REAR POUCH of biohazard bag and SEAL.
- 3 Place biohazard bag into SILVER INSULATED BAG > then pack into CSS/ ALCAT BOX. (Do not use FedEx Drop Box).
- 4 Place CSS/ ALCAT BOX INSIDE the FedEx Bag. In case of extreme summer temperature, use frozen ice packs.

Sample Collection Guide

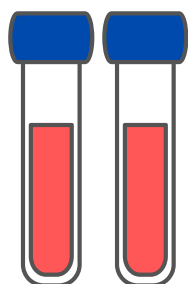

- ✓ For Full Array (Cellular Micronutrient Assay + REDOX/ APA = Cellular Nutrition Assay):
Blood Collection: WHOLE BLOOD
Requires immediate shipment.
2x Blue Top 3.8% Sodium Citrate Tube

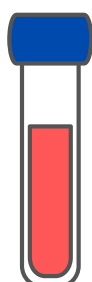

- ✓ For Individual Test Options (Cellular Micronutrient Assay ONLY or REDOX/ APA ONLY):
Same sample type as Full Array.
1x Blue Top 3.8% Sodium Citrate Tube

For Questions, Concerns, or Support: