

Step 1 - Preparation




- ✓ Ensure patient has fasted for 12 hours prior to collection.
- ✓ Fill out requisition form completely and legibly.
- ✓ Label vials: Full legal name, DOB, Date & Time of Collection.

Step 2 - Collect and Package Correctly




- 1 Place TUBES in SLEEVE > Place in BIOHAZARD BAG
- 2 Place REQUISITION FORM into REAR POUCH of biohazard bag and SEAL.
- 3 Ship out samples from Monday to Thursday ONLY. (US Clients)

Sample Collection Guide

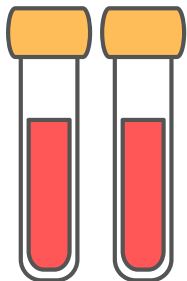

- ✓ For CardioPro Basic (Ship ASAP) | Serum Sample
2x Serum Separator Tube (SST)
 - Slowly invert SST Tubes 5-10 times.
 - Place the SST Tubes in a tube rack vertically to clot for no less than 30 minutes and no more than 60 minutes.
 - Centrifuge blood for 15 minutes at 3000-4000 rpm.
 - SST Tubes must have a balancing tube from it in the centrifuge with the same amount of liquid.

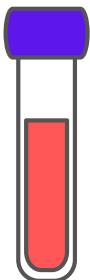

- ✓ For Hemoglobin A1C (Ship ASAP) | Whole Blood EDTA
1x Lavender Tube
 - Invert the lavender tube 8-10 times to allow the mixture of additives.

For Questions, Concerns, or Support:

 concierge@ownitcoaching.com

 +1 (510) 422-0209